

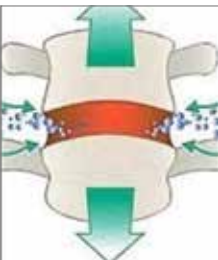
THERAMASTER

COMPONENTS



- ① Control Panel
- ② Cervical decompression
- ③ Cervical angle adjustment
- ④ Cervical one side
- ⑤ Height control
- ⑥ Emergency switch for patient
- ⑦ Up & down mobilization
- ⑧ Sacral mobilization
- ⑨ Pelvic mobilization
- ⑩ Lumbar angle adjustment
- ⑪ Lumbar decompression
- ⑫ Lumbar one side
- ⑬ Infrared / Up
- ⑭ Infrared / Down

THERAMASTER



CORE 1
• DISC DECOMPRESSION



CORE 2
• SPINE CORRECTION



CORE 3
• SPINE MOBILIZATION



CORE 4
• INTELLIGENT PROGRAM



www.spinemt.net

Manufactured by Shinhwa medical
THERAMASTER
MASTER OF SPINAL THERAPY

You can use it, depending on your purpose such as;

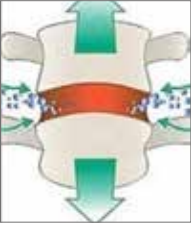
1. Basic chiropractic table + manipulation
(James M..Cox & James Cyriax Techniques)
2. Premium auto decompression with programmed Technics
3. Combination of 1 and 2 in the above.

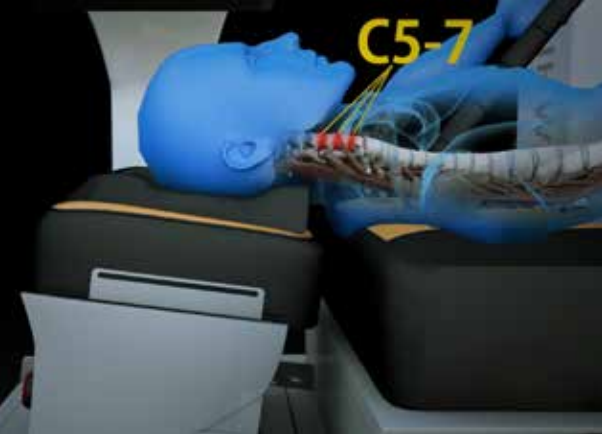
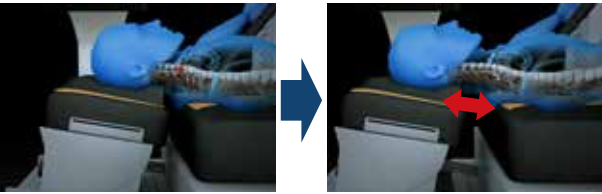
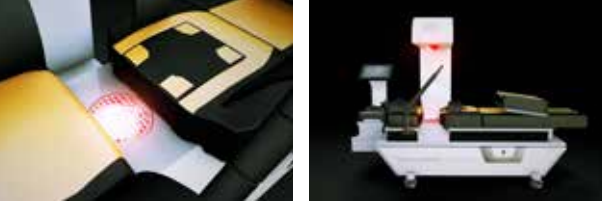


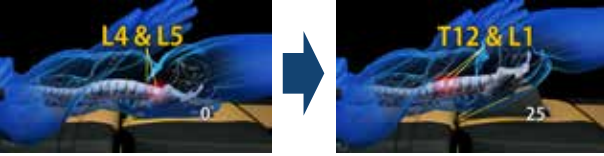
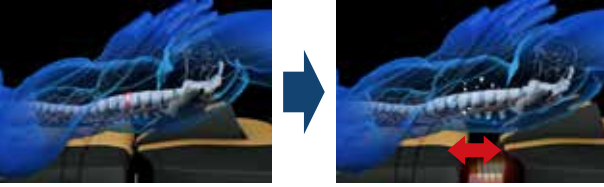
SHINHWA MEDICAL. Co., Ltd.
Address : #20, Cheongwon-gil, Ilgwang-myeon, Gijang-gun, Busan 46047 KOREA
Tel. +82-70-4946-4333 Fax. +82-51-724-6644
H.P. +82-10-2504-6655 E-mail. a8614493@naver.com www.spinemt.net



CORE THERAPY

	CORE 1 : DECOMPRESSION • Feedback System • Comfortable, Effective
	CORE 2 : CORRECTION • Target Set-up • Function : ①Unilateral ②Rotation ③Level
	CORE 3 : EXERCISE PRESCRIPTION • Intensity set-up : ①Low ②Medium ③High • Angle Set-up : 5 Steps
	CORE 4 : PROGRAM • Automatic Mode : Variously sophisticated Treatments
INDICATIONS	EFFECTS
• HNP & Degenerative Disc • Anti-colytic Spondylitis • Scoliosis & Stenosis • Joint Problems • Muscle Imbalance for rebuild	• Reduce pain • Increase vertebral flexibility • Reduce muscle spasm • Relieve nerve compression • Resorption of protruding discs

Part/Function	Image	Explanation
Cervical Target Angle		• Treatment principle (flexion-distraction) of Cox • Function to cure the cervical region to be treated through target angle adjustment • Adjustable from 0 ~ +15°
		• Functions according to medical examination and treatment principle • Functions to apply treatment according to the herniated direction of the disk • Up to 80mm adjustable for right & left
Cervical Decompression therapy		• Application of feedback system • Application of exercise function • Application of decompression function • Application of automatic program
Infrared Radiation therapy		• Up / Down Infrared ray is radiated to the rigid patient's body to relax the shortened muscles during the treatment. • The function maximizes the performance of the treatment. • Up infrared ray is optional.

Part/Function	Image	Explanation
Lumbar Target Angle		• Treatment principle (flexion-distraction) of Cox • Function to cure the lumbar region to be treated through target angle adjustment • Level adjustable from 0 ~ + 25°
Lumbar Unilateral application		• Functions according to medical examination and treatment principle • Functions to apply treatment according to the herniated direction of the disc • Up to 90mm adjustable for right & left
Lumbar Rotation application		• Treatment principle of Siriaks (orthopedic surgeon) • Treatment according to the type of escape of the disk • Adjust the angle of lumbar motion • Stop position from 0 to 18 ° for twisted spine & pelvis
Lumbar Decompression therapy		• Application of feedback system • Application of exercise function • Application of decompression function • Application of automatic program